


# LOW PRESSURE HOSE ASSEMBLIES EN ISO 5359

The low pressure hose assemblies, meeting the requirements of the Standard EN ISO 5359, are indicated for the use with the following medical gases:

- Oxygen (O<sub>2</sub>);
- Nitrous oxide (N<sub>2</sub>O);
- Medical air (AIR);
- Air for driving surgical tools (AIR-800);
- Vacuum (VAC).

These flexible hose assemblies can be used in a pressure range from 300 kPa to 1400 kPa for compressed medical gases and from 10 kPa to 100 kPa (absolute pressure) for vacuum. Here under some examples of suggested uses for the flexible hose assemblies:


- a.** Connection between terminal units and medical equipment (refer to Standards EN ISO 9170-1, EN ISO 80601-2-13 and EN ISO 80601-2-12);
- b.** Connection between a distribution pipeline system and one of its terminal units (refer to Standards EN ISO 7396-1 and EN ISO 11197);
- c.** Connection between two terminal units (refer to Standard EN ISO 7396-1);
- d.** Connection between an emergency supply source and a supply point for the emergency and the maintenance of a distribution pipeline system (refer to Standards EN ISO 10524-1 and EN ISO 7396-1);
- e.** Connection between an emergency supply point and a medical equipment (refer to Standards EN ISO 10524-1, EN ISO 80601-2-13 and EN ISO 80601-2-12).

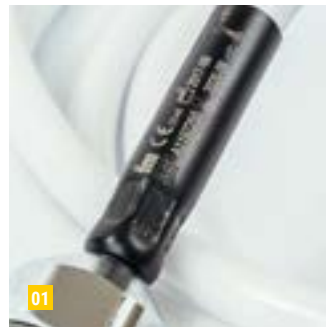

01



02

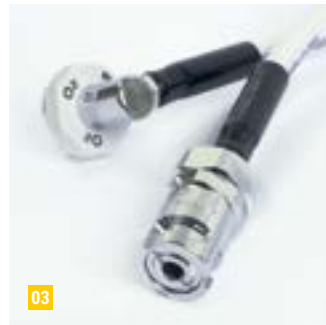

03

- 01 UNAMOVABLE CLAMPING SYSTEM
- 02 AVAILABLE FOR ALL MEDICAL GASES
- 03 DIFFERENT CONFIGURATIONS OF PROBES AND OUTLETS


**TECHNICAL SPECIFICATIONS | Low pressure hose assemblies EN ISO 5359**
**COMPRESSED MEDICAL GASES**
**HOSE SIZES**

Ø 11.5x5 mm

**VACUUM HOSE SIZES**

Ø 13.5x6.3 mm

**HOSE INNER MATERIAL**

black EPDM with antistatic treatment

**HOSE ARMOUR**

polyester

**HOSE COATING**

EPDM coloured according to the gas type

**COMPRESSED MEDICAL GASES**
**WORKING PRESSURES**

from 300 kPa to 1400 kPa

**VACUUM WORKING PRESSURES**

from 10 kPa to 100 kPa (absolute pressure)

**STANDARD HOSE LENGTHS**

1.5 m, 2.5 m, 4.5 m for gases

1.5 m, 2.5 m for vacuum