

EasyOX[®]

BUBBLING HUMIDIFIERS

The bubbling humidifiers for oxygen therapy series **EasyOX[®]** are devices used to increase the relative humidity in the treatments with medical oxygen, both for hospital and homecare applications.

Medical oxygen, normally supplied from hospital plants or from cylinders, is not provided with a sufficient degree of humidity to be physiologically tolerated by the patient's airways without collateral consequences. Hence the need to equip oxygen supply devices, such as the flowmeters, with suitable systems for the gas humidification during therapies.

The humidifiers for oxygen series **EasyOX[®]** are completely made of polycarbonate and they can be sterilized in autoclave at 121 °C for 15 min. The advanced design, the interchangeable snap threaded nut, allowing to fit the device in any situation, and the overall ergonomics make the humidifier **EasyOX[®]** particularly innovative. In conjunction with flowmeters **EasyFLOW[®]**, **Rs**, **Qmed[®]**, **EASYMED[®]** and **EasyMED[®] PLUS**, this humidifier guarantees the operators an easy use and great versatility, combined with an extreme economy.


The bottle hollow shape, besides being a design unique element, offers the operator an easy and safe handling and guarantees a higher humidification efficiency.



01



02

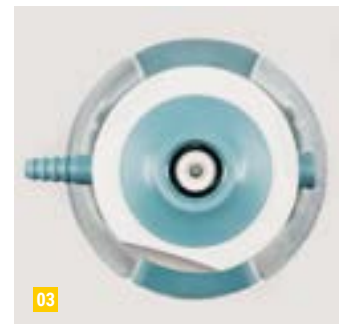

03



04

01 WITH SNAP INSERT, IDENTIFIED WITH COLOR-CODE, PROVIDED WITH DIFFERENT CONNECTIONS TO BE CHOSEN BY THE END USER

02 OPENING BY QUICK 1/12 TURN TO FACILITATE THE HUMIDIFIER FILLING AND CLEANING PROCEDURES

03 AN ERGONOMIC AND BIG SIZE ROTATING NUT ALLOWS TO EASILY CONNECT AND DISCONNECT THE HUMIDIFIER **EasyOX[®]** TO THE OXYGEN SUPPLY DEVICE. 360° ROLLING POSITIONING HOSE CONNECTOR Ø 6÷9 MM

04 DETAIL OF THE RELIEF VALVE

EXCELLENCE IN INNOVATION
BREAKTHROUGH TECHNOLOGY
INTELLIGENT DESIGN



flow-meterTM
ENGINEERING REVOLUTION IN MEDICAL TECHNOLOGY

REGISTERED OFFICE AND FACTORY ADDRESS: Flow Meter Spa - Via del Lino, 6 - 24040 Levate - Bergamo - ITALY
Ph. +39 035 594047 | Fax +39 035 59482 | info@flowmeter.it | www.flowmeter.it

VAT. N. IT01737040160


TECHNICAL SPECIFICATIONS | EasyOX®

INNER MAXIMUM VOLUME	190 ml
SIZE (LxWxH)	89x84x161 mm
WEIGHT	0.13 Kg
MAXIMUM APPLICABLE PRESSURE VALUE	80 kPa
MAXIMUM APPLICABLE FLOW VALUE	10 L/min.
INLET CONNECTION	9/16" UNF EN 13544-2 F. • 1/4" ISO 3253 F. M 12x1.25 F. • 3/8" ISO 3253 F.
OUTLET CONNECTION	hose connector Ø 6÷9 mm